

GLOVES FOR ISOLATOR

ถุงมือยาวสำหรับตู้เตรียมตัวอย่าง



PIERCAN drybox gloves

Reference Guide:
Standard Models



To find the reference number
of the glove that meets your needs,
follow these steps.



- Medical
- Electronic
- Laboratory Testing
- Food Science Research
- Battery Vacuum Glove Box
- Pharmaceutical Research Labs
- Nuclear
- Humidity Control
- Industrial Research



บริษัท 3เอ อินสตรูเมนต์ จำกัด

เลขที่ 200/61 หมู่ 5 โทรทัศน์ 0 2046 0691-2

ตำบลบางรักน้อย 0 2046 2608

อำเภอเมืองนนทบุรี โทรสาร 0 2045 2608

จังหวัดนนทบุรี 11000 www.3ainstrument.com

มือถือ./Mobile Phone. +668 1508 1168, +666 3045 5536, +669 4539 5052

3A INSTRUMENT CO.,LTD.

200/61 Moo.5

Tambon Bangruknoi,

Amphur Muang Nonthaburi,

Nonthaburi 11000 Thailand

Tel. +66 2046 0691-2

+66 2046 2608

Fax. +66 2045 2608

E-mail: Sales@3ainstrument.com

Contact@3ainstrument.com



Bechtop Glove Box : SUS80 Series



3A Instrument Co., Ltd.

for control environment with close system.

These basic compact Dry glove boxes provide low-humidity, low-oxygen Free controlled environment when used with Close system . Includes inlet-outlet gas port for inert gas condition testing . **Glove boxes are ideal choice for laboratory testing cost system in such industries as Laboratory, Research Testing, pharmaceutical manufacturing, Cosmetic, Semiconductor, Batteries, Electronic , Solar cell and University research.** Dry Glove Box system designed for easy transfer of samples and close system for dry glove box.


Features:

- **304 Stainless steel thickness 2.0 mm. Tilt - up front window clear Acrylic thickness 12 mm with large opening. Safety with shock-up**
- can be cleaning inside system with the inert gas port.
- Two glove port diameter 200 mm.
- Dry glove boxes provide low-humidity, controlled environment when used with Close system.
- Glove Box design for negative pressure and positive pressure testing.
- Pass box allow convenient sample parts transfer.



Application

- Close system testing for Laboratories of pharmaceuticals, Biological, Chemical.
- Weighing samples of moisture sensitive drug powders.
- Experiments that require low-humidity conditions.
- Operation of various instruments requiring environmental control, Equipment Enclosure such as Automatic titration, Volumetric and Culometric Karl Fischer Titration.
- **Anaerobic Test Glove Box.**
- **Negative or positive pressure glove box.**
- **Vacuum Glove box.**
- **Biosafety Glove Box.**
- **Benchtop Glove Box for control Environment Equipment Enclosure.**
- **Gas control system nitrogen purge for low –humidity.**
- **Temperature and Dehumidity control with PLC Controller.**
- **Basic Research Glove Box.**

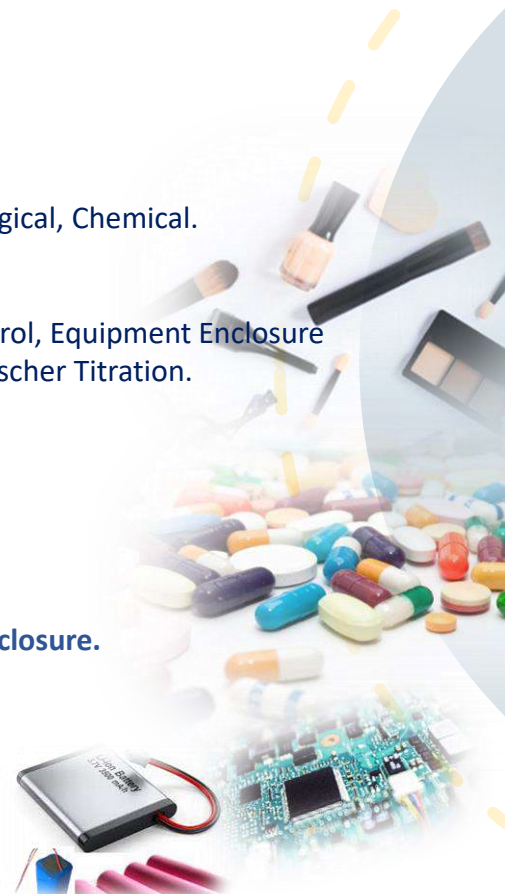

Basic Glove Box

- Dry glove boxes provide low-humidity, controlled environment when used with Close system.
- Glove Box design for negative pressure and positive pressure testing.
- Environment Equipment Enclosure for application Microscope , Balance.



Chamber

- Easy open and operate system .
- Glove Port ID 200 mm for use Long gloves Diameter 220 mm, Length 800 mm.
- Gloves Box Model S80PLUS With ... is special system. There are 2 systems, namely the system that filters clean air (Can choose filters to the needs of the work) and automatic systems that control the temperature and humidity inside the gloves box. .
- Gas (N₂, O₂, or Air) in port Push-to Connect fitting and gas-out bleeding valve. Flow meter for gas control system. Glove port with O-Ring 1 Set (2 pcs)



Pass Box

- The AIR BOX is Accessories transfer samples into the glove box Size W350 x D350 x H350 mm. Side open window. The 300 x 300 mm. opening allows for parts to be easily placed inside the glove box with door slide.
- Gas (N₂, O₂, or Air) in port Push-to Connect fitting and gas-out bleeding valve. Flow meter for gas control system. Glove port with O-Ring 1 Set (2 pcs)




Specifications:

Model	SUS80 -30	SUS80 -31	SUS80 -3123	SUS80 -3124	SUS80 -3125	SUS80 -3127	SUS80 -3128	SUS80 -3129
Material Stainless Steel SUS304.	Stainless Steel SUS304							
Material SUS304 Thickness :	Thickness 2 mm							
External Dimensions (Main Chamber)	W 800 x D 600 x H 600 mm.							
Lighting 15 W .	Yes.							
Gas Port : for inert gas such as nitrogen, argon, air, etc.	Yes. (2 pcs.)							
Analog Pressure Gauge : 0-10 bar	Yes.							
Door open big side : form Acrylic Material Thickness 12 mm .	W 500 x D 720 mm.							
Glove Port Diamention (2 pcs) :	ID 200 mm							
Pass Box for sample transfer : - Size : W350 x D350 x H350 mm - Inlet/Outlet Gas Port , 2 pcs - Analog Pressure 0-10 bar, 1 pcs	-	Yes.	Yes.	Yes.	Yes.	Yes.	Yes.	Yes./ 2 set
Carbon Filter with Blower	-	-	Yes.	-	-	-	-	
Hepa Filter with Blower	-	-	-	Yes.	-	-	-	
Auto Dry for low humidity	-	-	-	-	Yes.	-	-	
Temperature & Dehumidity Control - Thermoelectric Cooler - Anti-Condensation Heater - Dehumidifier	-	-	-	-	-	Yes	-	-
Programmable process Controller - Control Temperature step - Dehumidity Control	-	-	-	-	-	-	Yes 99 step Yes	-
PLC Controller	-	-	-	-	-	-	-	Yes.
O ₂ Purification units.	Options.							
Power Socket : 220V,50-60 Hz,	Yes. (2 Channel.)							
Gloves: Port ID 220 mm,Length 800 mm	Yes. (Options) Material : Natural Rubber /Neoprene/CSM/BHP/PUR/EPDM							





Drawing size dimensions

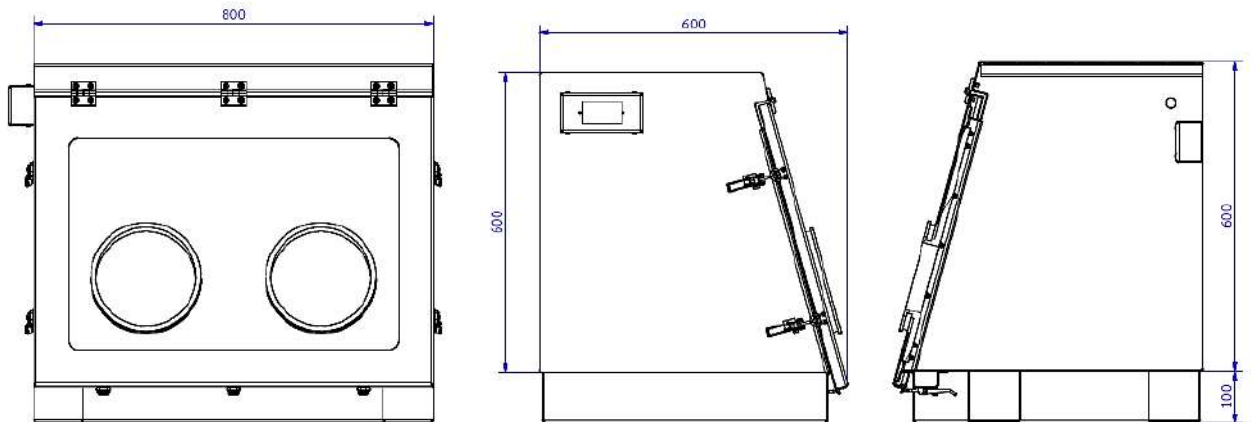

Main Chamber


- **Door working area :**
Clear Acrylic : thickness 12 mm.
- **External Chamber Diamention :**
W 800 x D 600 x H 600 mm.
- **Main Chamber : Stainless Steel ,**
Thickness 2 mm
- **Open Door internal**
Diamention : 700 x 450 mm.

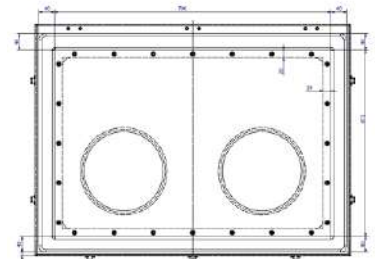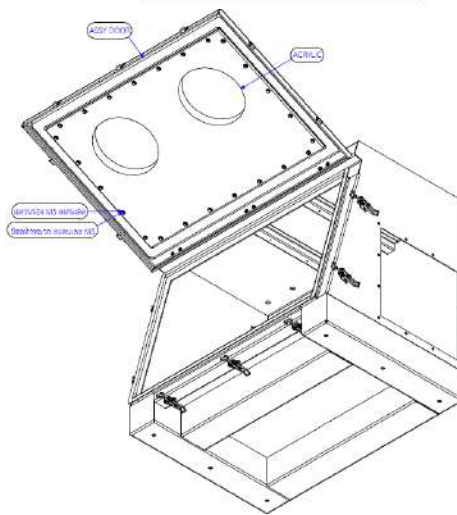

Pass Box

- **Door Clear Acrylic :**
thickness 25 mm.
- **External Diamention :**
350 x 350 x 350 mm.
- **Chamber Stainless steel 304**
- **Open Door internal**
Diamention : 300 x 300 mm.
- **Inlet/Outlet for inert gas**
Port 2 pcs.

Door Open Transfer sample

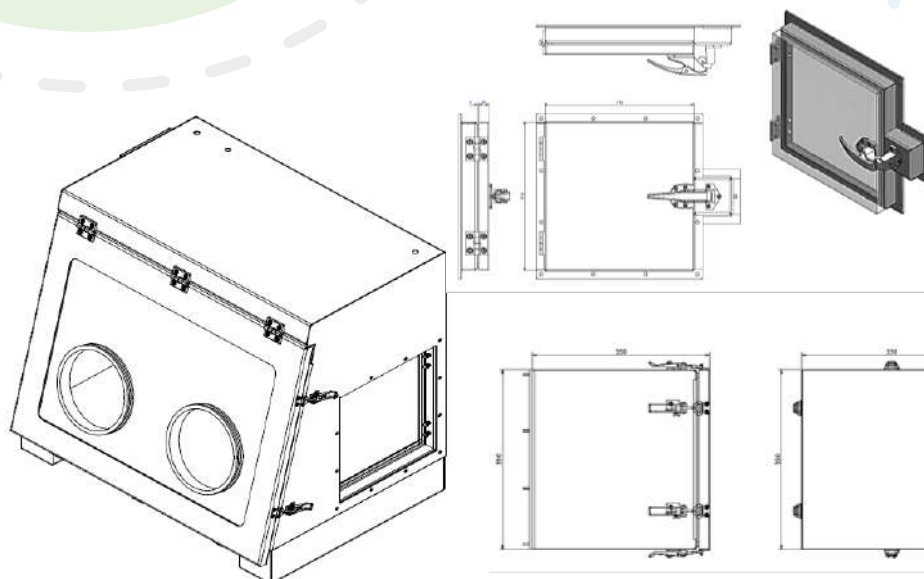

Glove Box : SUS100 Series

for control environment with close system.

These basic compact Dry glove boxes provide low-humidity, low-oxygen Free controlled environment when used with Close system . Includes inlet-outlet gas port for inert gas condition testing . **Glove boxes are ideal choice for laboratory testing cost system in such industries as Laboratory , Research Testing , pharmaceutical manufacturing , Cosmetic , Semiconductor , Batteries , Electronic , Solar cell and University research.** Dry Glove Box system designed for easy transfer of samples and close system for dry glove box.


Features:

- **304 Stainless steel thickness 2.0 mm. Tilt-up front window clear Acrylic thickness 12 mm with large opening. Safety with shock-up**
- can be cleaning inside system with the inert gas port.
- Two glove port diameter 200 mm.
- Dry glove boxes provide low-humidity, controlled environment when used with Close system.
- Glove Box design for negative pressure and positive pressure testing.
- Pass box allow convenient sample parts transfer.



Application

- Close system testing for Laboratories of pharmaceuticals, Biological, Chemical.
- Weighing samples of moisture sensitive drug powders.
- Experiments that require low-humidity conditions.
- Operation of various instruments requiring environmental control, Equipment Enclosure such as Automatic titration, Volumetric and Coulometric Karl Fischer Titration.
- **Anaerobic Test Glove Box.**
- **Negative or positive pressure glove box.**
- **Vacuum Glove box.**
- **Biosafety Glove Box.**
- **Benchtop Glove Box for control Environment Equipment Enclosure.**
- **Gas control system nitrogen purge for low –humidity.**
- **Temperature and Dehumidity control with PLC Controller.**
- **Basic Research Glove Box.**

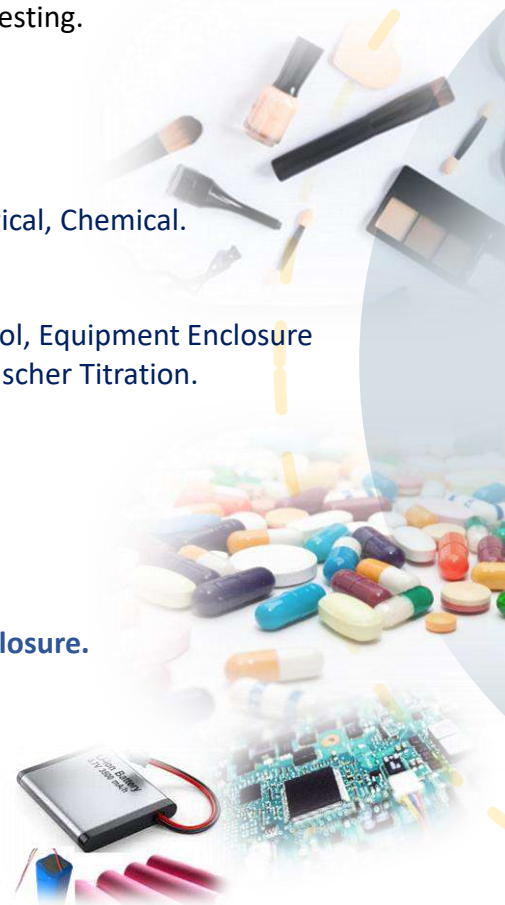

Basic Glove Box

- Dry glove boxes provide low-humidity, controlled environment when used with Close system.
- Glove Box design for negative pressure and positive pressure testing.
- Environment Equipment Enclosure for application Microscope , Balance.



Chamber

- Easy open and operate system .
- Glove Port ID 200 mm for use Long gloves Diameter 220 mm, Length 800 mm.
- Gloves Box Model S80PLUS With ... is special system. There are 2 systems, namely the system that filters clean air (Can choose filters to the needs of the work) and automatic systems that control the temperature and humidity inside the gloves box. .
- Gas (N₂, O₂, or Air) in port Push-to Connect fitting and gas-out bleeding valve. Flow meter for gas control system. Glove port with O-Ring 1 Set (2 pcs)



Pass Box

- The AIR BOX is Accessories transfer samples into the glove box Size W350 x D350 x H350 mm. Side open window. The 300 x 300 mm. opening allows for parts to be easily placed inside the glove box with door slide.
- Gas (N₂, O₂, or Air) in port Push-to Connect fitting and gas-out bleeding valve. Flow meter for gas control system. Glove port with O-Ring 1 Set (2 pcs)




Specifications:

Model	SUS100 -30	SUS100 -31	SUS100 -3123	SUS100 -3124	SUS100 -3125	SUS100 -3127	SUS100 -3128	SUS100 -3129
Material Stainless Steel SUS304.	Stainless Steel SUS304							
Material SUS304 Thickness :	Thickness 2 mm							
External Dimensions (Main Chamber)	W 1000 x D 700 x H 800 mm.							
Working Area Dimensions	W 1000 x D 700 x H 700 mm.							
Lighting 15 W .	Yes.							
Gas Port : for inert gas such as nitrogen, argon, air, etc.	Yes. (2 pcs.)							
Analog Pressure Gauge : 0-10 bar	Yes.							
Door open big side : form Acrylic Material Thickness 12 mm .	W 600 x D 900 mm.							
Glove Port Diamention (2 pcs) :	ID 200 mm							
Pass Box for sample transfer : - Size : W350 x D350 x H350 mm - Inlet/Outlet Gas Port , 2 pcs - Analog Pressure 0-10 bar, 1 pcs	-	Yes.	Yes.	Yes.	Yes.	Yes.	Yes.	Yes./ 2 set
Carbon Filter with Blower	-	-	Yes.	-	-	-	-	
Hepa Filter with Blower	-	-	-	Yes.	-	-	-	
Auto Dry for low humidity	-	-	-	-	Yes.	-	-	
Temperature & Dehumidity Control - Thermoelectric Cooler - Anti-Condensation Heater - Dehumidifier	-	-	-	-	-	Yes	-	-
Programmable process Controller - Control Temperature step - Dehumidity Control	-	-	-	-	-	-	Yes 99 step Yes	-
PLC Controller	-	-	-	-	-	-	-	Yes.
O2 Purification units.	Options.							
Power Socket : 220V,50-60 Hz,	Yes. (2 Channel.)							
Gloves: Port ID 220 mm,Length 800 mm	Yes. (Options) Material : Natural Rubber /Neoprene/CSM/BHP/PUR/EPDM							




Drawing dimensions

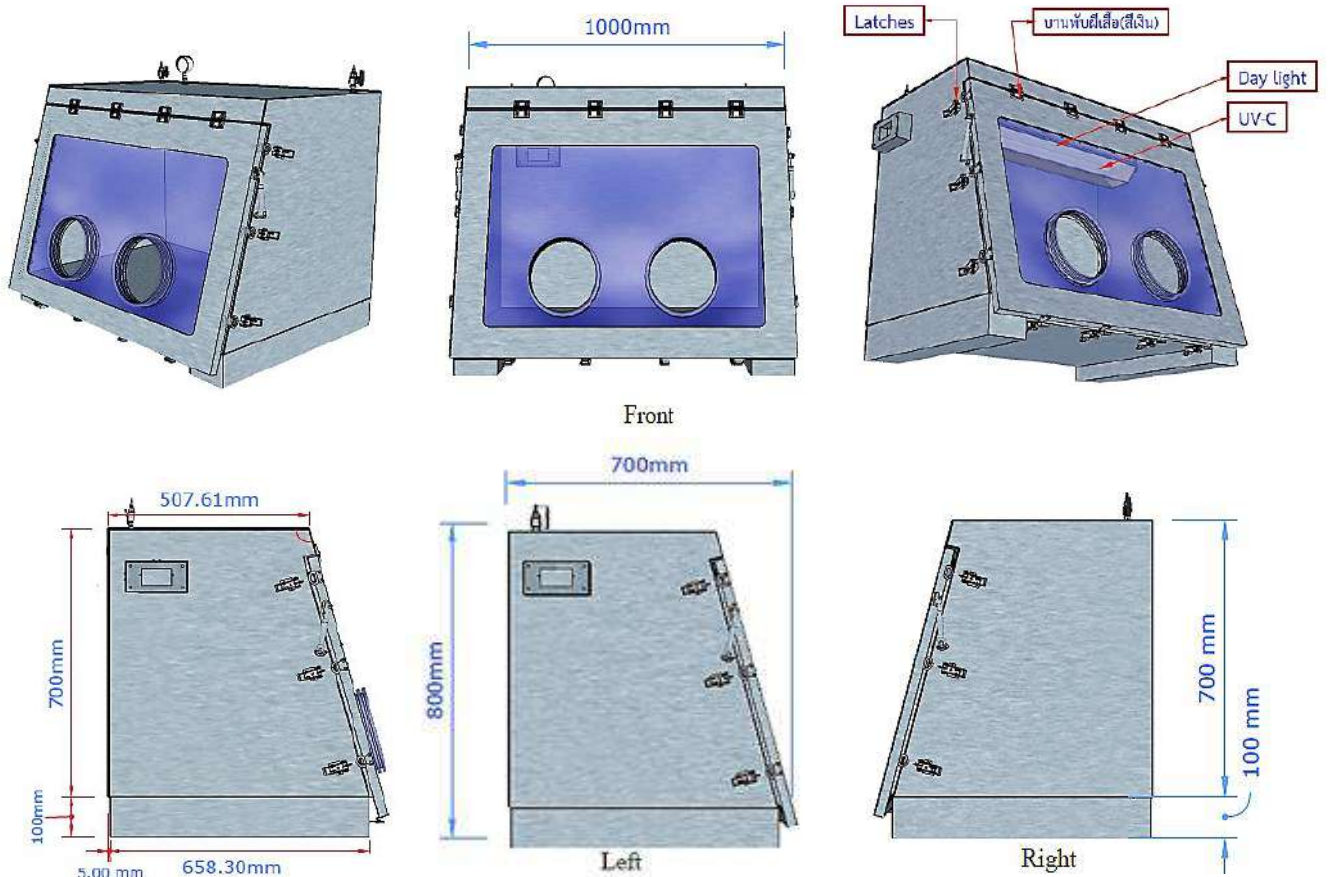






Main Chamber

Easy Maintenance and Testing

- **External Chamber Dimension :**
W 1000 x D 700 x H 800 mm.
- **Working Area Dimension :**
W 1000 x D 700 x H 700 mm.
- **Main Chamber : Stainless Steel Thickness : 2 mm**
- **Open Door internal Dimension : 900 x 600 mm.**

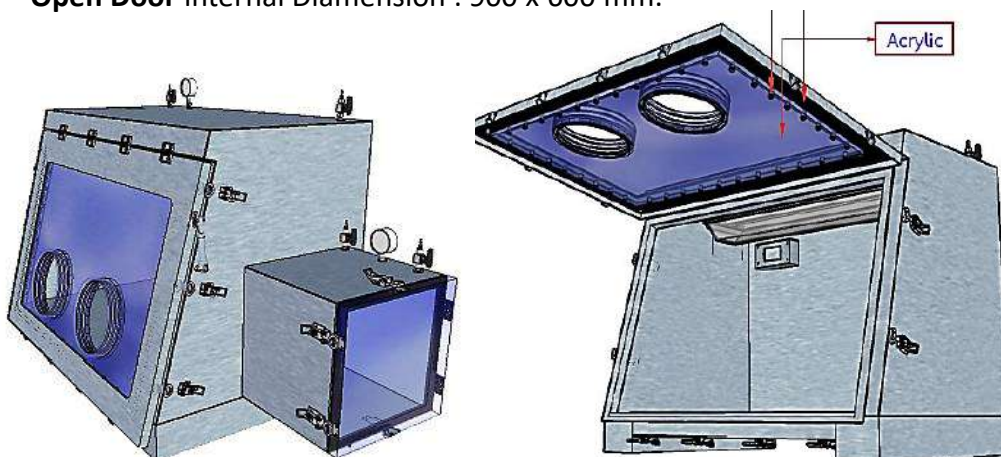


Drawing dimensions

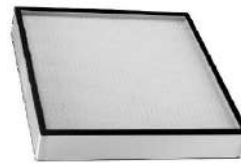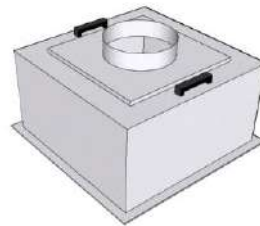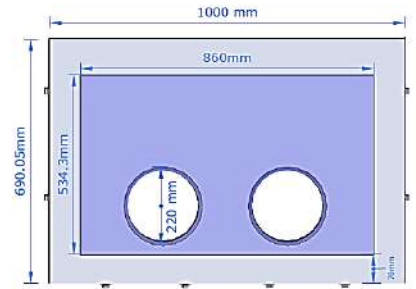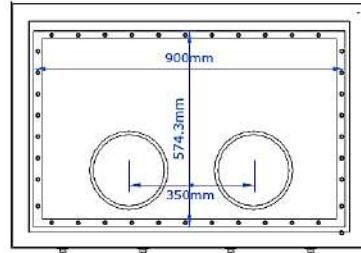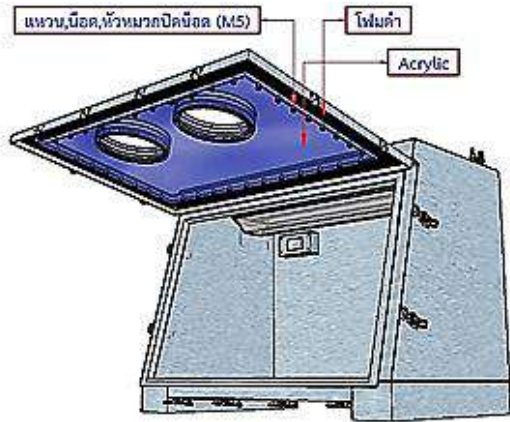

Hepa Filter






- **HEPA-Filtered with blower Glove Boxes** have inlet and outlet HEPA filters and provide barrier to protect the operator from hazardous airborne particulates and powders.
- **Carbon- filter with Blower Glove Box** Helps prevent dust and chemicals.

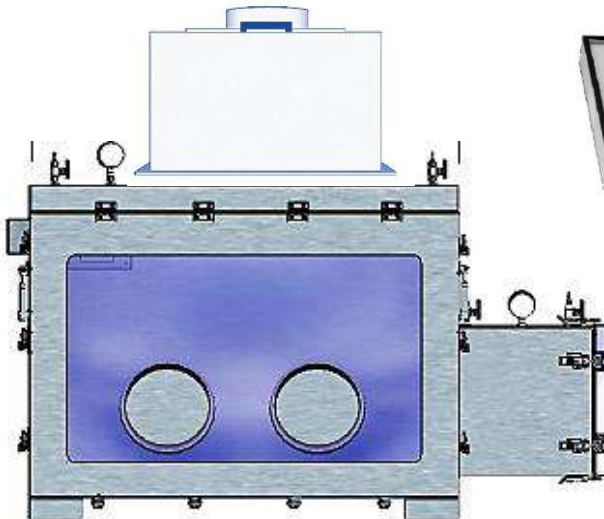

Pass Box

- **Door Clear Acrylic :** thickness 25 mm.
- **External Diamention :** 350 x 350 x 350 mm.
- **Chamber Stainless steel 304**
- **Open Door internal** Diamention : 300 x 300 mm.
- **Inlet/Outlet for inert gas** Port 2 pcs.

Pass box for Transfer sample

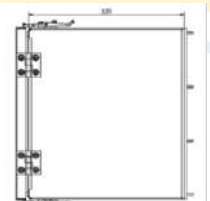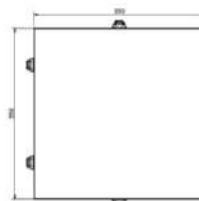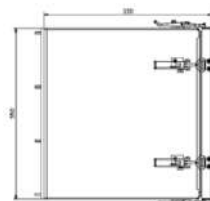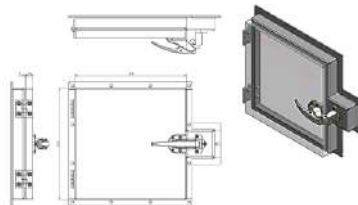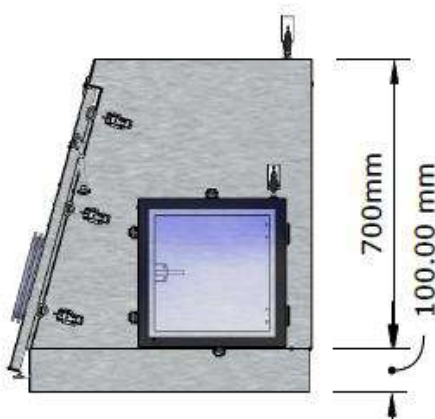

Glove Box : SUS120 Series

for control environment with close system.

These basic compact Dry glove boxes provide low-humidity, low-oxygen Free controlled environment when used with Close system . Includes inlet-outlet gas port for inert gas condition testing . ***Glove boxes are ideal choice for laboratory testing cost system in such industries as Laboratory , Research Testing , pharmaceutical , Cosmetic , Semiconductor , Batteries , Electronic , Solar cell and University research.*** Dry Glove Box system designed for easy transfer of samples and close system for dry glove box.


Features:

- **304 Stainless steel thickness 2.0 mm. Tilt-up front window clear Acrylic thickness 12 mm with large opening. Safety with shock-up**
- can be cleaning inside system with the inert gas port.
- Two glove port diameter 200 mm.
- Dry glove boxes provide low-humidity, controlled environment when used with Close system.
- Glove Box design for negative pressure and positive pressure testing.
- Pass box allow convenient sample parts transfer.



Application

- Close system testing for Laboratories of pharmaceuticals, Biological, Chemical.
- Weighing samples of moisture sensitive drug powders.
- Experiments that require low-humidity conditions.
- Operation of various instruments requiring environmental control, Equipment Enclosure such as Automatic titration, Volumetric and Coulometric Karl Fischer Titration.
- **Anaerobic Test Glove Box.**
- **Negative or positive pressure glove box.**
- **Vacuum Glove box.**
- **Biosafety Glove Box.**
- **Benchtop Glove Box for control Environment Equipment Enclosure.**
- **Gas control system nitrogen purge for low –humidity.**
- **Temperature and Humidity control.**
- **Programmable process Controller Temperature and Humidity ..**
- **Basic Research Glove Box.**

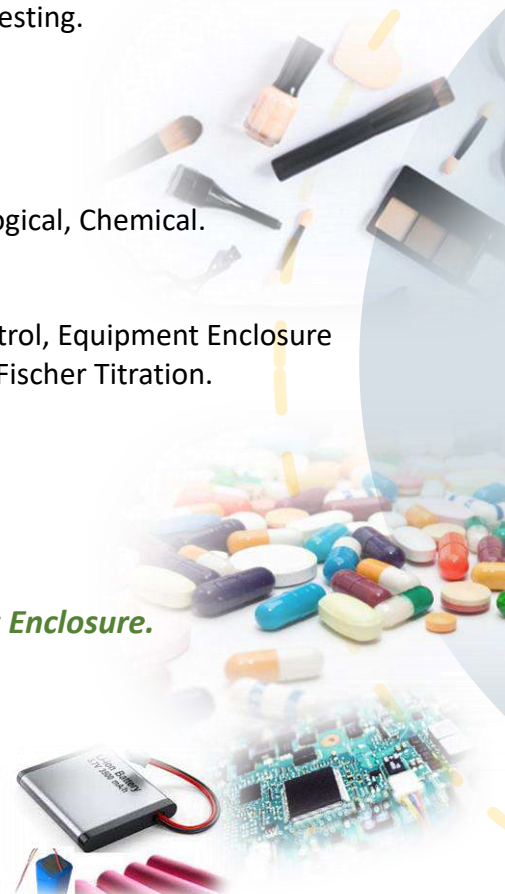

Basic Glove Box

- Dry glove boxes provide low-humidity, controlled environment when used with Close system.
- Glove Box design for negative pressure and positive pressure testing.
- Environment Equipment Enclosure for application Microscope , Balance.



Chamber

- Easy open and operate system .
- Glove Port ID 200 mm for use Long gloves Diameter 220 mm, Length 800 mm.
- Gloves Box Model S80PLUS With ... is special system. There are 2 systems, namely the system that filters clean air (Can choose filters to the needs of the work) and automatic systems that control the temperature and humidity inside the gloves box. .
- Gas (N₂, O₂, or Air) in port Push-to Connect fitting and gas-out bleeding valve. Flow meter for gas control system. Glove port with O-Ring 1 Set (2 pcs)



Pass Box

- The AIR BOX is Accessories transfer samples into the glove box Size W350 x D350 x H350 mm. Side open window. The 300 x 300 mm. opening allows for parts to be easily placed inside the glove box with door slide.
- Gas (N₂, O₂, or Air) in port Push-to Connect fitting and gas-out bleeding valve. Flow meter for gas control system. Glove port with O-Ring 1 Set (2 pcs)




Specifications:

Model	SUS120 -30	SUS120 -31	SUS120 -3123	SUS120 -3124	SUS120 -3125	SUS120 -3127	SUS120 -3128	SUS120 -3129
Material Stainless Steel SUS304.	Stainless Steel SUS304							
Material SUS304 Thickness :	Thickness 2 mm							
External Dimensions (Main Chamber)	W 1200 x D 800 x H 900 mm.							
Working Area Dimensions	W 1200 x D 800 x H 800 mm.							
Lighting 15 W .	Yes.							
Gas Port : for inert gas such as nitrogen, argon, air, etc.	Yes. (2 pcs.)							
Analog Pressure Gauge : 0-10 bar	Yes.							
Door open big side : form Acrylic Material Thickness 12 mm .	W 1100 x D 900 mm.							
Glove Port Diamention (2 pcs) :	ID 200 mm							
Pass Box for sample transfer : - Size : W350 x D350 x H350 mm - Inlet/Outlet Gas Port , 2 pcs - Analog Pressure 0-10 bar, 1 pcs	-	Yes.	Yes.	Yes.	Yes.	Yes.	Yes.	Yes./ 2 set
Carbon Filter with Blower	-	-	Yes.	-	-	-	-	
Hepa Filter with Blower	-	-	-	Yes.	-	-	-	
Auto Dry for low humidity	-	-	-	-	Yes.	-	-	
Temperature & Dehumidity Control - Thermoelectric Cooler - Anti-Condensation Heater - Dehumidifier	-	-	-	-	-	Yes	-	-
Programmable process Controller - Control Temperature step - Dehumidity Control	-	-	-	-	-	-	Yes 99 step Yes	-
PLC Controller	-	-	-	-	-	-	-	Yes.
O2 Purification units.	Options.							
Power Socket : 220V,50-60 Hz,	Yes. (2 Channel.)							
Gloves: Port ID 220 mm,Length 800 mm	Yes. (Options) Material : Natural Rubber /Neoprene/CSM/BHP/PUR/EPDM							






Drawing dimensions

Easy Maintenance and Testing


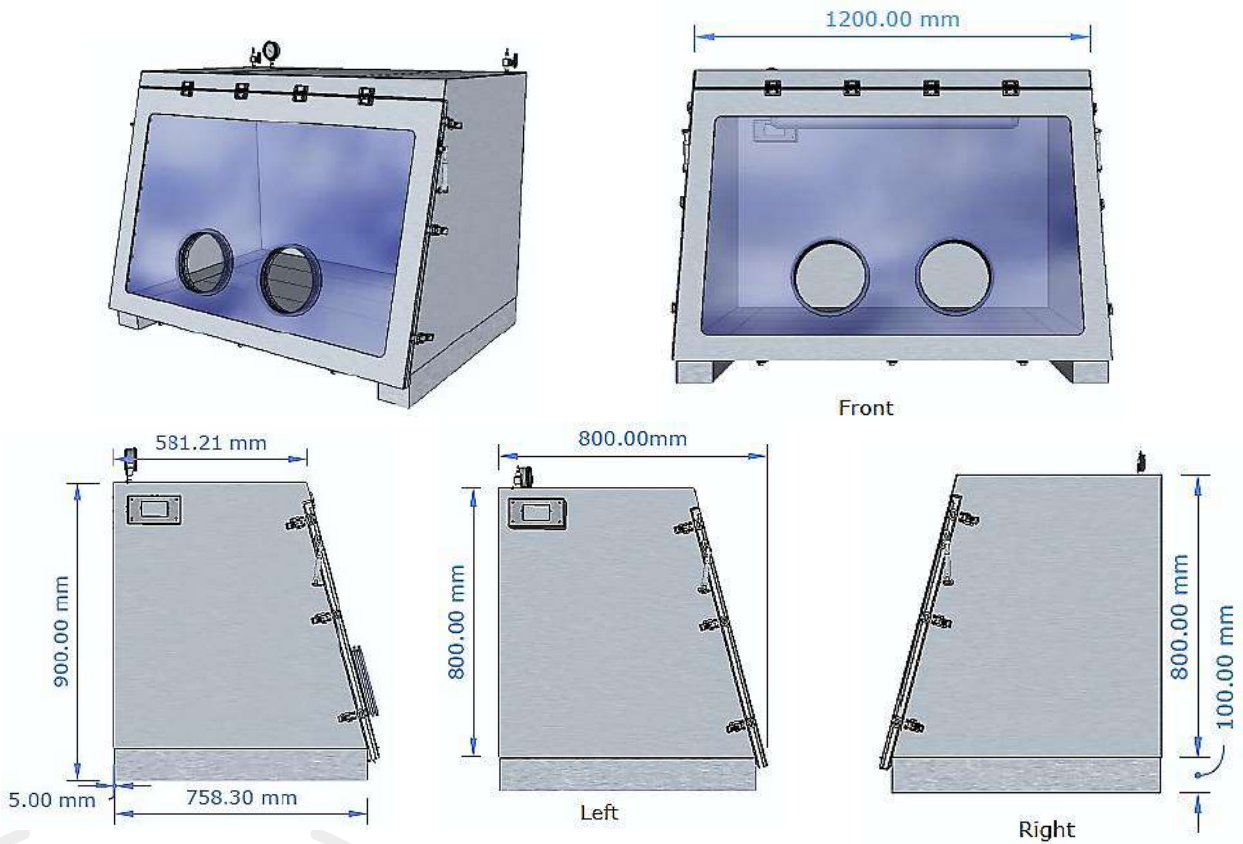


Main Chamber

- External Chamber Dimension :
W 1200 x D 800 x H 900 mm.
- Working Area Dimension :
W 1200 x D 800 x H 800 mm.
- Main Chamber : Stainless Steel
Thickness : 2 mm
- Open Door internal Dimension :
1100 x 700 mm.

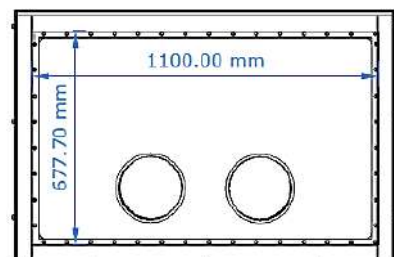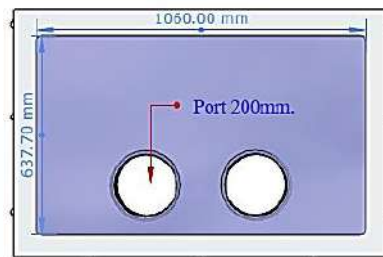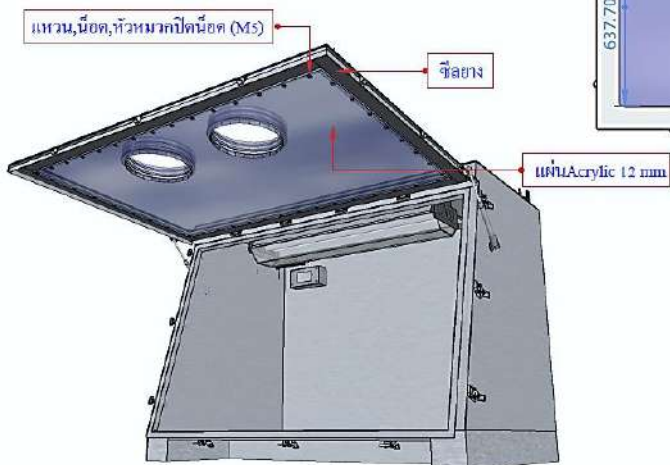

Drawing dimensions

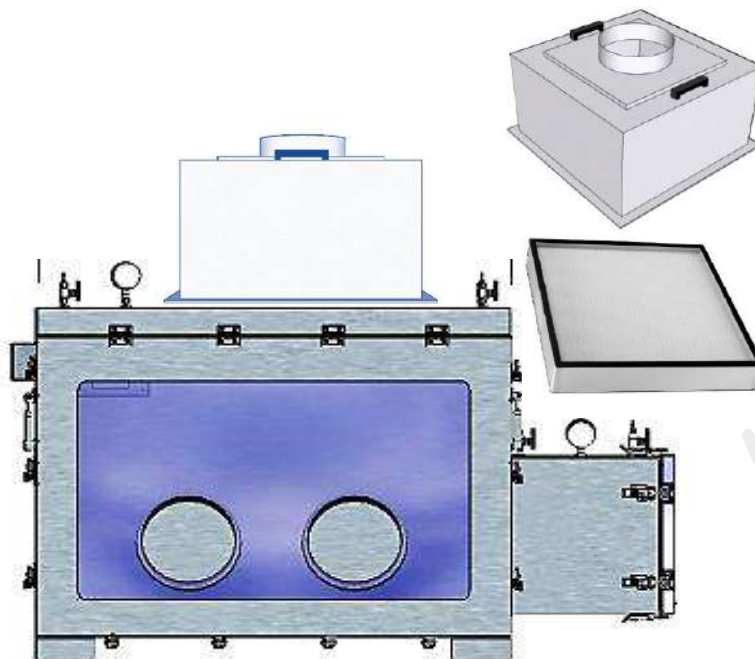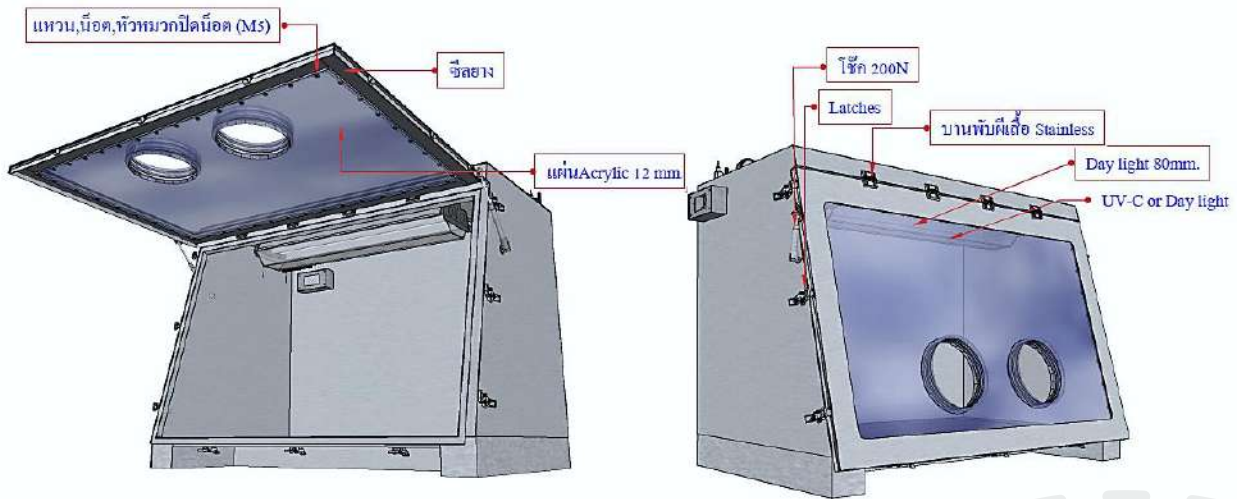


Hepa Filter

- **HEPA-Filtered with blower Glove Boxes** have inlet and outlet HEPA filters and provide barrier to protect the operator from hazardous airborne particulates and powders.
- **Carbon- filter with Blower Glove Box** Helps prevent dust and chemicals.

Pass Box

- **Door Clear Acrylic** : thickness 25 mm.
- **External Diamention** : 350 x 350 x 350 mm.
- **Chamber Stainless steel 304**
- **Open Door internal Diamention** : 300 x 300 mm.
- **Inlet/Outlet for inert gas** Port 2 pcs.

Door Open Transfer sample

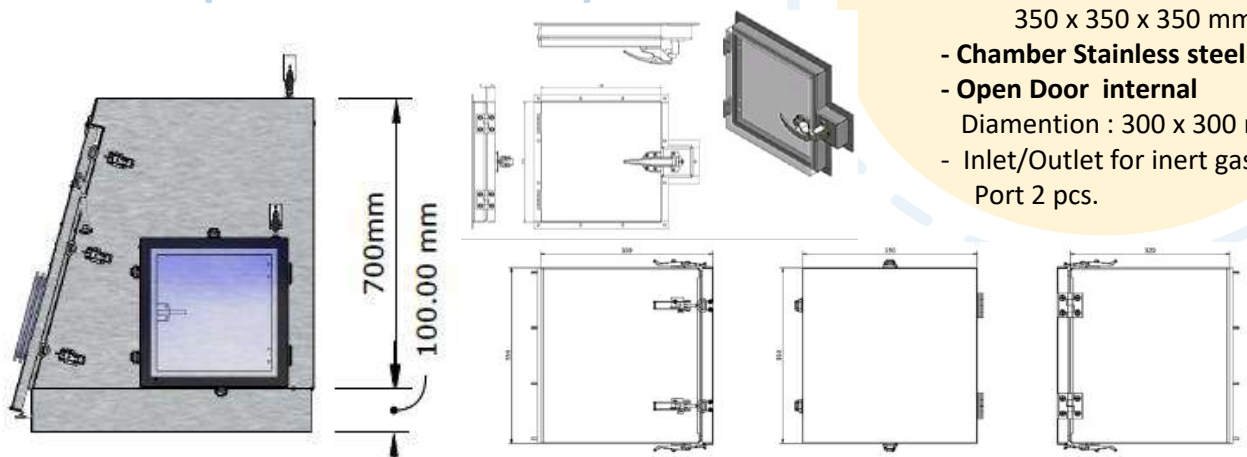



Molecular sieve



Humidity Temperature Datalogger



POM Glove Port



Gloves for ISOLATOR PIERCAN FRANCE



Oil Free Laboratory Chemical Resistant Vacuum




EDPM O-ring for Glove Port



Air to air thermoelectric Cooler



Vacuubrand MZ 2C NT 1.2CFM Chemistry Diaphragm Pump 110V - Chemical Resistance

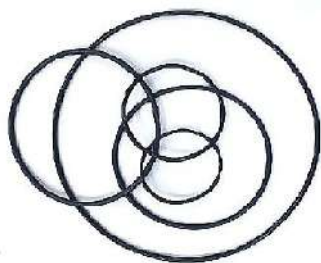

Digital Temperature Controller

Technical Dry Gloves for Glove Box

Technical Material gloves for all industrial and laboratory testing.



PIERCAN Glove is high quality product form France. The range of gloves that PIERCAN offers is very wide, comprising dry box gloves, sleeve, and isolator glove.



CSM

Chlorosulfonated Polyethylene glove (CSM).

- Resistance to ozone, UV and Sterilizing agents (VHP*).
- Excellent mechanical properties.
- Resistance to ionizing radiation and sterilizing process (Gamma and Beta Radiation).

VHP* (Vaporous hydrogen peroxide)



Black EPDM

Ethylene Propylene Diene Monomer glove.

- Antistatic
- Excellent Flexibility and dexterity.
- Excellent chemical properties, good mechanical properties.
- Composition in accordance with the FDA positive list (CFR 21 § 177. 2600).
- Resistance to AUTOCLAVE, GAMMA & VHP Sterilization 75 cycles of 30 minutes at 121 °C.



Natural Rubber

Natural rubber glove (polyisoprene)

- Excellent flexibility and dexterity
- Good resistance to alcohols, reducing acids and diluted bases.

DRY GLOVE

for the Pharmaceutical and Cosmetics manufacturing.



Neoprene

Polychloroprene glove

- Good mechanical properties.
- Good general resistance to chemical products.
- Very good resistance to sterilisation agents (hydrogen peroxide) and radiation sterilisation (Gamma and Beta rays).
- Good technical and cost-effective compromise
- Self- extinguishing.



POLYURETHANE

POLYURETHANE glove

- Excellent mechanical properties (punctures, tears, abrasion)
- Highly resistant to ozone and UV rays



POLYURETHANE / CSM

Two layer polyurethane and chlorosulfonated polyethylene glove.

- A glove that combines mechanical properties (polyurethane) with chemical resistance (CSM)
- Good resistance to chemical sterilising agents (CSM side)
- High resistance to ozone and UV. Good resistance to ionising radiation.
- Excellent resistance to disinfectants.



High Butyl Properties

Polyisoprene isobutylene glove.

- Electrostatic dissipative glove in accordance with EN 16350.
- Highly impermeable to liquids and gases
- Good ability to withstand chemical products
- Highly resistant to ozone and UV rays
- Good flexibility and dexterity
- Antistatic (compliant with European standard EN 16350-2014)



บริษัท 3เอ อินสตรูเมนต์ จำกัด

ผู้ผลิตตู้เตรียมตัวอย่าง (Glove Box) และจัดจำหน่ายถุงมือสำหรับตู้ ISOLATOR
นำเข้าจากประเทศฝรั่งเศส

ติดต่อสอบถาม

คุณปิยะนุช โพธิ์หล้า	Piyanuch@3ainstrument.com	08 1508 1168
คุณประทุมพร คุณล้าน	Sales@3ainstrument.com	09 4559 5052
คุณอภิสิทธิ์ ละครลำ	Apisit@3ainstrument.com	06 1152 9338
คุณวัชรพงษ์ สิทธิน้อย	Info@3ainstrument.com	09 3543 5103



LINE ID : @3ainstrument